

(19)



Europäisches Patentamt
European Patent Office
Office européen des brevets



(11) Publication number:

0 361 468 B1

(12)

EUROPEAN PATENT SPECIFICATION

- (45) Date of publication of patent specification: **26.05.93** (51) Int. Cl.⁵: **G01N 33/92, G01N 33/68, G01N 33/577, G01N 33/538**
- (21) Application number: **89117937.6**
- (22) Date of filing: **28.09.89**

(54) **A kit for the detection of denatured lipoproteins.**

(30) Priority: **29.09.88 JP 244987/88**
06.09.89 JP 231893/89

(43) Date of publication of application:
04.04.90 Bulletin 90/14

(45) Publication of the grant of the patent:
26.05.93 Bulletin 93/21

(84) Designated Contracting States:
DE FR GB

(56) References cited:
EP-A- 0 257 778

**CHEMICAL ABSTRACTS, vol. 108, no. 13,
March 28, 1988, Columbus, OH (US); T.SAITO
et al., p. 325, no. 109122***

**CHEMICAL ABSTRACTS, vol. 107, no. 17,
October 26, 1987, Columbus, OH (US);
S.KSAHARA et al., p. 360, no. 150520m***

**CHEMICAL ABSTRACTS, vol. 106, no. 26,
June 29, 1987, Columbus, OH (US); T.MORI et
al., p. 375, no. 219579m***

(73) Proprietor: **JAPAN IMMUNO RESEARCH LAB-
ORATORIES CO., LTD.**
17-5, Sakae-cho
Takasaki-shi Gunma-ken 370(JP)

(72) Inventor: **Adachi, Masakazu**
3493-9, Ishihara-machi
Takasaki-shi Gunma-ken(JP)
Inventor: **Saitou, Toshimitu**
1928-340, Oaza Munataka Gunma-cho
Gunma-gun Gunma-ken(JP)
Inventor: **Murata, Kimiyo**
B-101, 239-7, Fujizuka-machi
Takasaki-shi Gunma-ken(JP)
Inventor: **Tamura, Aki**
2-9-6, Ohtone-cho
Maebashi-shi Gunma-ken(JP)

(74) Representative: **Wächtershäuser, Günter, Dr.**
Tal 29
W-8000 München 2 (DE)

Note: Within nine months from the publication of the mention of the grant of the European patent, any person may give notice to the European Patent Office of opposition to the European patent granted. Notice of opposition shall be filed in a written reasoned statement. It shall not be deemed to have been filed until the opposition fee has been paid (Art. 99(1) European patent convention).

EP 0 361 468 B1

Description

BACKGROUND OF THE INVENTION5 Field of the Invention:

The present invention relates to a kit and a method for detecting denatured lipoproteins, and, more particularly, to a kit and a method for detecting denatured lipoproteins useful for the clinical investigation and diagnosis of abnormal lipid metabolism which can be a cause of heart diseases (e.g. cardiac
 10 infarction), diabetes, and various types of atherosclerosis diseases such as cerebrovascular disorders (e.g. cerebral hemorrhage).

Description of the Background Art:

15 Abnormal lipid metabolism can be defined as abnormality in either synthesis, secretion, or degradation of lipoproteins. Abnormal lipoprotein metabolism includes congenital abnormality and acquired abnormality.

Abnormal lipid metabolism is considered to be a cause of various diseases such as heart diseases (e.g. cardiac infarction) and various types of atherosclerosis diseases such as cerebrovascular disorders (e.g. cerebral hemorrhage), and to have significant influence on the progress and aftercare of these diseases.
 20 Detection of abnormal lipid metabolism is thus very important for the clinical investigation and diagnosis of these diseases.

A method most widely accepted for detecting abnormal lipid metabolism is that using the total cholesterol value in blood as a parameter. With this method, the total cholesterol value in blood is measured by a conventional method (e.g. an enzymatic method), and abnormal lipid metabolism is judged to be
 25 present if the value falls outside the normal range which was established based on clinical and experiential knowledge. In another method used in recent years lipoproteins in blood, i.e., combined lipid - apoproteins, are fractionated by specific gravities, and a cholesterol value for high specific gravity lipoproteins (high density lipoproteins; HDL) is measured. The HDL cholesterol value, or the value obtained by subtracting the HDL cholesterol value from the total cholesterol value (the atherosclerosis index), is used as a parameter for
 30 detecting abnormal lipid metabolism.

According to the results of many immunological studies and clinical knowledge accumulated up to the present time, however, these conventional methods of detecting abnormal lipid metabolism do not always exactly reflect the incidence of abnormal lipid metabolism and have only a limited clinical significance. This is because in these methods the total or HDL cholesterol values, which are not necessarily associated with
 35 abnormal lipid metabolism, are used. Although the high value of total cholesterol value can be one of the risk factors of ischemic heart diseases [e.g. *Med. Clin. North Am.*, 58, 363-379 (1974)], it does not directly relate to the diseases. Occurrence of the diseases in groups with normal lipid metabolism or healthy conditions in groups with highly abnormal lipid metabolism according to these detection methods is not rare. Detection using as a parameter a HDL cholesterol value or using an atherosclerosis index based on
 40 the HDL cholesterol value does not always accurately reflect the incidence of disease *Arteriosclerosis* 14, (4), 931-936 (1986)].

Many reports have recently surfaced discussing the significance on actual abnormal lipid metabolism of denatured lipoproteins which are produced when lipoproteins such as very low density lipoproteins (VLDL) undergo certain types of actions *in vivo* [e.g. Hui, D. Y., Innerarity, T. L. and Mahley, R. W., *J. Bio. Chem.*,
 45 259, 860-869 (1984); Gonen, et al., *Diabetes*, 30, 875 (1981)].

Detection of such denatured lipoproteins can therefore more accurately detect the incidence of abnormal lipid metabolism and is thus considered to be clinically very useful as a diagnostic means for determining diseases involving abnormal lipid metabolism, especially for the early discovery of the diseases and for monitoring the progress or curing effects of the diseases. Details of such denatured lipoproteins are,
 50 however, not yet elucidated at the present time and there are no known methods for their detection.

The present inventors have conducted extensive studies in order to develop a method and a kit for accurately detecting abnormal lipid metabolism which can detect various diseases involving abnormal lipid metabolism and can monitor the progress or curing effects of such diseases. As a result, the present inventors found that a portion obtained by excluding lipoproteins constituting apoprotein A-I (Apo A-I -
 55 containing lipoprotein) and lipoproteins constituting apoprotein B-100 (Apo B-100 - containing lipoprotein) from total plasma lipoproteins corresponded to the above denatured lipoproteins. The inventors further found that the Apo A-I-containing lipoprotein and the Apo B-100-containing lipoprotein could be removed when a blood sample was incubated with affinity gels containing specific antibodies, thus ensuring

an accurate and ready measurement or detection of the amount of denatured lipoproteins in the sample. Such a finding has led to the completion of the present invention.

SUMMARY OF THE INVENTION

5

Accordingly, an object of the present invention is to provide a kit for detecting denatured lipoproteins which comprises as an essential component a mixture of an affinity gel with which a monoclonal antibody for apoprotein A-I is coupled and an affinity gel with which a monoclonal antibody for apoprotein B-100 is coupled.

10

Another object of the present invention is to provide a method for measuring an amount of denatured lipoproteins which comprises: adding to a sample a mixture of an affinity gel with which a monoclonal antibody for apoprotein A-I is coupled and an affinity gel with which a monoclonal antibody for apoprotein B-100 is coupled; removing apoprotein A-I and apoprotein B-100 in the sample by combining the apoprotein A-I and apoprotein B-100 with the affinity gels; and measuring the amount of lipoproteins in

15

the residue. Other objects, features and advantages of the invention will hereinafter become more readily apparent from the following description.

BRIEF DESCRIPTION OF THE DRAWINGS

20

Fig. 1 is a graph showing the correlation between the amounts of denatured lipoproteins measured using the kit of the present invention and the plasma cholesterol values.

Fig. 2 is a graph showing the plasma cholesterol values of a hyperlipemia patient before and after treatment.

25

Fig. 3 is a graph showing the amounts of denatured lipoprotein of a hyperlipemia patient before and after treatment.

DETAILED DESCRIPTION OF THE INVENTION AND PREFERRED EMBODIMENTS

30

The kit of the present invention is designed to measure the amount of denatured lipoproteins, i.e., lipoproteins excluding Apo A-I-containing lipoprotein and Apo B-100-containing lipoprotein, from the total lipoproteins in blood. According to the present invention it is possible to detect abnormal lipid metabolism and to clinically investigate and diagnose various diseases involving abnormal lipid metabolism by using the kit. Since the above denatured lipoproteins are rarely present in blood when lipid metabolism

35

is normal and increase depending on the degree of abnormality of lipid metabolism, their amount in blood accurately reflects any abnormality of the lipid metabolism. Thus, an accurate measurement and detection of abnormal lipid metabolism can be ensured by using the kit of the present invention.

The kit of the present invention as well as the measurement of denatured lipoproteins or the detection of abnormal lipid metabolism by using the kit are described in detail in the following.

40

A sample used in the detection by using the kit of the present invention is blood. Serum or plasma collected while the subject is under fasting conditions is especially preferable as samples. Samples are prepared for the test according to conventional methods after having been collected from subjects.

45

Since lipoproteins exist in blood as a combination of a lipid such as a cholesterol, triglyceride, phospholipid, free fatty acid, or the like and an apoprotein, either the total amount of lipids or the total amount of apoproteins can be used for the measurement of the amount of lipoproteins. The cholesterol values for the lipids or for the apoproteins can be determined by an enzymatic method or the like.

50

A mixture of an affinity gel with which a monoclonal antibody for apoprotein A-I (anti-apoprotein A-I antibody) is coupled and an affinity gel with which a monoclonal antibody for apoprotein B-100 (anti-apoprotein B-100 antibody) is coupled is used as is or packed in a tube which can be shaken or centrifuged. Alternatively, it is desirable to equilibrate the gels in advance with a suitable buffer, for example, with a 0.1 M Tris-HCl (pH 7.4) + 0.15 M NaCl buffer. A conventional preservative such as sodium azide or the like can be added to the buffer.

55

An affinity gel used in the present invention can be prepared by using a conventional method. Specifically, the anti-apoprotein A-I antibody and the anti-apoprotein B-100 antibody are activated by cyanogen bromide and each is coupled with an affinity gel. The gels are mixed and equilibrated with a suitable buffer solution. Typical examples of antibodies which can be used are anti-human apo A-I monoclonal antibody and anti-human apo B-100 monoclonal antibody (both produced by Japan Immunoresearch Lab. Co.). Antibodies are by no means limited to these antibodies in this invention. A typical

example of an affinity gel is Sepharose 4B (trade name, product of Pharmacia Co.) activated by cyanogen bromide or the like in advance. A coupling reaction can be carried out according to a method conventionally used for coupling reactions of the above type of gels and proteins, usually at pH range of 8–10 using an antibody in an amount of about 1–10 mg for 1 ml of a gel at room temperature (20–25 °C) within 2 hours.

There are no specific limitations as to the proportion of anti-apoprotein A-I antibody and anti-apoprotein B-100 antibody in the above gel mixture which is the essential component of the kit of this invention. Usually, an equivalent or more, preferably 1–2 times by weight, of anti-apoprotein B-100 antibody is used for anti-apoprotein A-I antibody. It is desirable to adjust the pH value of the gel mixture to about 7–7.5.

If required, a stabilizer such as sucrose, a bovine serum protein, or the like, or a preservative can be added to the kit of the present invention containing the gel mixture as an essential component. A preservative is selected from compounds which will not adversely affect the result obtained by the use of the kit. A typical example is diluted sodium azide. Glycerol or a water-soluble glycerol derivative, alcohol, glycol, glycol ether, or the like can also be added to the kit.

According to a preferred embodiment for detecting denatured lipoproteins by using the kit of the present invention, the affinity gel mixture which is the essential component of the kit is placed in a suitable container or test tube and equilibrated with a suitable buffer. A prescribed amount of serum to be tested is added to the container or the tube and allowed to react at room temperature by being left to stand still. A cholesterol coloring reagent is added to the supernatant. After incubation for a prescribed period of time, an absorbance of the supernatant is measured to determine the amount of lipids or apoproteins contained therein. The amount of denatured lipoproteins in the serum sample is then determined. Usually, about 1/4–1/16 by weight, preferably about 1/8 by weight, of a sample serum is used for the gel mixture of the present invention.

Denatured lipoproteins can be readily measured by the use of the kit of the present invention. The amount of the denatured lipoproteins thus measured can be used as a parameter for detecting abnormal lipid metabolism. Thus, the kit of the present invention provides a ready and effective means for the detection of abnormal lipid metabolism. Since the amount of denatured lipoproteins in blood accurately reflects the abnormality of lipid metabolism and its degree, the kit of the present invention is very useful for the discovery of the above-mentioned diseases and for the monitoring of the progress or curing effects of such diseases. It is especially useful for the accurate and effective detection of heart diseases such as cardiac infarction and various types of obesities which have not previously been detectable by conventional detection methods.

Other features of the invention will become apparent in the course of the following description of the exemplary embodiments which are given for illustration of the invention.

EXAMPLES

Example 1

(1) Blood was collected with the addition of EDTA from a patient while fasting and centrifuged (2,500–3,000 rpm) to obtain plasma. The triglyceride content of the plasma was measured by an enzymatic method using TG-555 (trade name, product of Kyowa Medics Co., Ltd.). The cholesterol content was also determined by an enzymatic method using Mercko Test CHO (trade name, product of Kanto Chemical Co., Ltd.) in a 2 ml/10 μ l amount of plasma. The mixture was incubated for 10 minutes at 37 °C. The absorbance was measured at 365 nm.

(2) The above plasma (20 μ l) was charged into a test tube and to this was added 600 μ l of a reaction fluid containing an anti-apo A-I monoclonal antibody and an anti-apo B-100 monoclonal antibody, both coupled with an affinity gel [Sepharose 4B activated by cyanogen bromide; antibodies (as proteins) in 10 mg units were coupled with 1 ml of the gel and each gel was mixed with an equivalent amount of 0.01 M Tris-HCl buffer (pH 7.4)]. The tube was shaken for 30 minutes and then allowed to stand for 15 minutes. 2 ml of Mercko Test CHO color reagent (trade name, product of Kanto Chemical Co., Ltd.) was added to the supernatant and the mixture was incubated for 20 minutes at 37 °C. The amount of cholesterol was determined from the absorbance measured at 365 nm.

(3) As a result, triglycerides in plasma were found to have a higher correlation than cholesterol. From this, the possibility was confirmed of detecting denatured lipoproteins reflecting a chylomicron remnant which is an indicator of exogenous abnormal lipid metabolism and abnormal VLDL (β -) which is an indicator of endogenous abnormal lipid metabolism.

(4) The plasma cholesterol values and the denatured lipoprotein values of healthy persons were determined using the same procedures as above. Fig. 1 shows their correlations, in which the amounts of plasma cholesterol (mg/dl) are plotted along the ordinate and the amounts of denatured lipoproteins are plotted along the abscissa (mg/dl).

Fig. 1 demonstrates that in healthy subjects the denatured lipoprotein values are constant having values in the close vicinity of 0 irrespective of the cholesterol values.

(5) For comparison, the denatured lipoprotein values and cholesterol values are determined by the following method without using the above method of the present invention.

In this method, a column packed with an affinity gel coupled with anti-apo A-I monoclonal antibody and a column packed with an affinity gel coupled with anti-apo B-100 monoclonal antibody were used to determine each cholesterol value for apoprotein A-1-containing lipoproteins (absorbed fractions) and apoprotein B-100-containing lipoproteins (absorbed fractions) in plasma. The cholesterol values for denatured lipoproteins were then determined by subtracting the sum of these cholesterol values from the total plasma cholesterol value.

Specifically, each 250 μ l of plasma from the same healthy subjects as in (4) above was charged for absorption into the columns packed with an affinity gel coupled with anti-apo A-I monoclonal antibody and an affinity gel coupled with anti-apo B-100 monoclonal antibody. The columns were washed with 20 ml of a 0.1 M PBS buffer (pH 7.2) and eluted using 10 ml of a mixed solution of 1 M acetic acid and 0.5 M NaCl. A 1/2-fold 3.5 M Tris-HCl buffer was added to the eluant. To a 1 ml aliquot of this mixture was added 2 ml of Mereko Test CHO (trade name, product of Kanto Chemical Co., Ltd.). After incubating at 37°C for 30 minutes, the cholesterol value was measured by absorption at 365 nm.

A parameter (P) was calculated from the cholesterol values measured as above according to the following formula:

$$\text{Parameter (\%)} = 1 - \frac{B + C}{A} \times 100 \quad (I)$$

wherein A is the total plasma cholesterol value (mg/dl), B is the cholesterol value (mg/dl) of apoprotein A-1-containing lipoproteins separated from the plasma, and C is the cholesterol value (mg/dl) of apoprotein B-100-containing lipoproteins separated from the plasma.

The results are shown in the columns designated as "Comparison Measurement" in Table 1, in which the results obtained in (4) above by the use of the kit of the present invention are also presented (in the columns designated as "Invented Measurement").

Table 1

Subject No.	Plasma Cholesterol Value		Denatured Lipoprotein Value	
	Invented Measurement	Comparison Measurement	Invented Measurement	Comparison Measurement (P)
1	169	159	0	-16 (-10)
2	119	133	0	11 (8)
3	169	188	0	28 (15)
4	166	157	0	13 (8)
5	159	140	0	8 (6)
6	146	219	0	2 (1)

The results shown in the above table show that the measurement using the kit of the present invention ensures a more accurate determination of denatured lipoprotein values with smaller fluctuations than the conventional method. The invented method is simple and easy to use and provides excellent results.

Example 2

500 mg (2 tablets each containing 250 mg) of Probucol (trade name, produced by Dow Pharmaceutical Co.) was given twice each day in the morning and evening for 98 days to a subject who had a high triglyceride value (TG > 150) and a high total cholesterol value (TC > 250) in the determination in Example 1 (1) and was thus judged to be hyperlipidemic. The triglyceride value and the total plasma cholesterol value as well as the denatured lipoprotein value according to the method of the present invention of the subject were measured in the same manner as in Example 1 before and after the 98-day period.

The triglyceride value changed from 425 (the value before treatment) to 389 (the value after treatment). Changes in the plasma cholesterol value and the denatured lipoprotein value are respectively shown in Figs. 2 (1) and (2). These data demonstrate that the kit and the method of the present invention ensure accurate determinations of denatured lipoproteins corresponding to symptoms and are useful for the diagnosis of diseases.

Example 3

The total cholesterol values (T-CHO), triglyceride values (TG), HDL cholesterol values (HDL-C), and denatured lipoprotein values (DLP-CHO) were measured on each 5-member group suffering from diabetes, ischemic heart disease, and cerebrovascular disorder. T-CHO, TG, and HDL-C were measured by the same methods as described in Example 1. DLP-CHO was measured by the heparin-manganese coupling precipitation method using HDL-cholesterol Wako (trade name, product of Wako Pure Chemical Co., Ltd.). The results are shown in Table 2.

TABLE 2

Subject No.	T-CHO (mg/dl)	TG (mg/dl)	HDL-C (mg/dl)	DLP-CHO (mg/dl)
Diabetes				
1	333	202	45	6
2	241	202	43	14
3	167	247	51	5
4	192	527	40	36
5	189	179	36	8
Ischemic heart disease				
1	241	187	47	12
2	296	812	33	81
3	198	127	42	8
4	143	94	59	9
5	331	190	43	12
Cerebrovascular disorder				
1	198	118	33	10
2	248	484	37	31
3	244	138	48	5
4	178	199	37	6
5	249	54	40	4

Example 4

Denatured lipoprotein values in the plasma of the hyperlipidemic patients were measured according to the method described in Example 1. The patients were grouped into three types as defined by WHO. The results are shown in Table 3.

TABLE 3

Types of hyperlipidemic patients	Number of patients	Denatured lipoprotein values *
Ila	31	2.1 ± 3.9 (mg/dl CHO)
IIb	53	20.1 ± 25.0
IV	16	7.7 ± 6.7

* Mean value ± S.D.

Claims

1. A kit for detecting denatured lipoproteins which comprises as an essential component a mixture of an affinity gel with which a monoclonal antibody for apoprotein A-I is coupled and an affinity gel with which a monoclonal antibody for apoprotein B-100 is coupled.
2. A method for measuring the amount of denatured lipoproteins which comprises: adding to a sample a mixture of an affinity gel with which a monoclonal antibody for apoprotein A-I is coupled and an affinity gel with which a monoclonal antibody for apoprotein B-100 is coupled; removing apoprotein A-I and apoprotein B-100 in the sample by combining the apoprotein A-I and apoprotein B-100 with the affinity gels; and measuring an amount of lipoproteins in the residue.

Patentansprüche

1. Testsatz zum Bestimmen denaturierter Lipoproteine, umfassend als wesentlichen Bestandteil ein Gemisch eines Affinitätsgels, an das ein monoklonaler Antikörper für Apoprotein A-I gekuppelt ist und ein Affinitätsgel, an das ein monoklonaler Antikörper für Apoprotein B-100 gekuppelt ist.
2. Verfahren zum Messen einer Menge denaturierten Lipoproteins, umfassend: Zugabe eines Gemisches aus einem Affinitätsgel, an das ein monoklonaler Antikörper für Apoprotein A-I gekuppelt ist und einem Affinitätsgel, an das ein monoklonaler Antikörper für Apoprotein B-100 gekuppelt ist zu einer Probe, Entfernen des Apoprotein A-I und Apoprotein B-100 der Probe durch Kombinieren des Apoproteins A-I und Apoprotein B-100 mit den Affinitätsgelen; und Messen der Lipoproteinmenge im Rest.

Revendications

1. Kit de détection de lipoprotéines dénaturées, qui comprend, comme composant essentiel, un mélange d'un gel d'affinité avec lequel un anticorps monoclonal pour l'apoprotéine A-I est couplé, et un gel d'affinité avec lequel un anticorps monoclonal pour l'apoprotéine B-100 est couplé.
2. Procédé de mesure de la quantité de lipoprotéines dénaturées, qui comprend :
 - l'addition, à un échantillon, d'un mélange d'un gel d'affinité avec lequel un anticorps monoclonal pour l'apoprotéine A-I est couplé, et d'un gel d'affinité avec lequel un anticorps monoclonal pour l'apoprotéine B-100 est couplé ;
 - l'élimination de l'apoprotéine A-I et de l'apoprotéine B-100 dans l'échantillon par combinaison de l'apoprotéine A-I et de l'apoprotéine B-100 avec les gels d'affinité ; et
 - la mesure de la quantité de lipoprotéines dans le résidu.

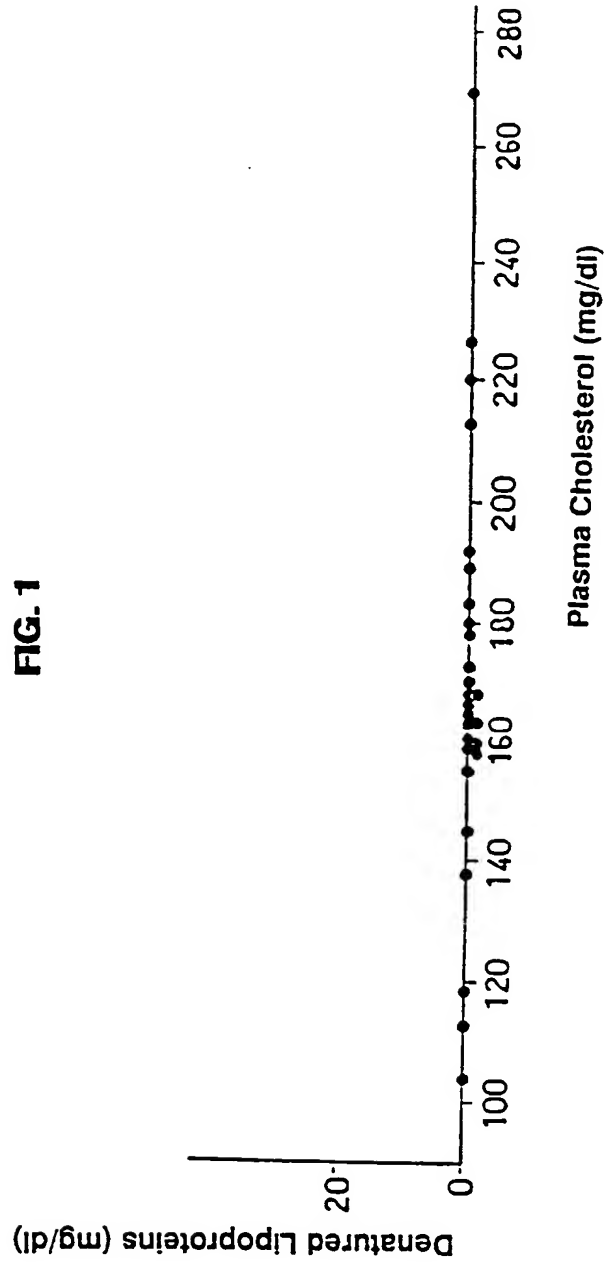


FIG. 2

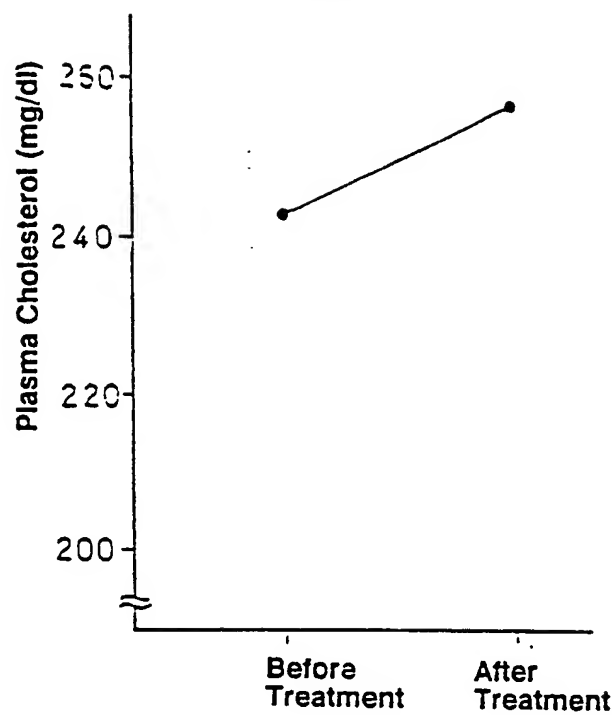


FIG. 3

